

Sample name: **Pyridoxine HCl**  
 Assay name: **UV-metric pKa**  
 Assay ID: **18E-22004**  
 Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

Experiment start time: **5/22/2018 10:11:10 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Results

pKa 1 **4.83**  
 pKa 2 **8.87**  
 RMSD **0.007 0.003 0.007**  
 Chi squared **0.0453**  
 PCA calculated number of pKas **2**  
 Average ionic strength **0.155 M**  
 Average temperature **24.9°C**  
 Analyte concentration range **61.7 µM to 57.8 µM**

Number of pKas source **Predicted**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.466 to 12.552**


## Warnings and errors

Errors None  
 Warnings None

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			









## Graphs

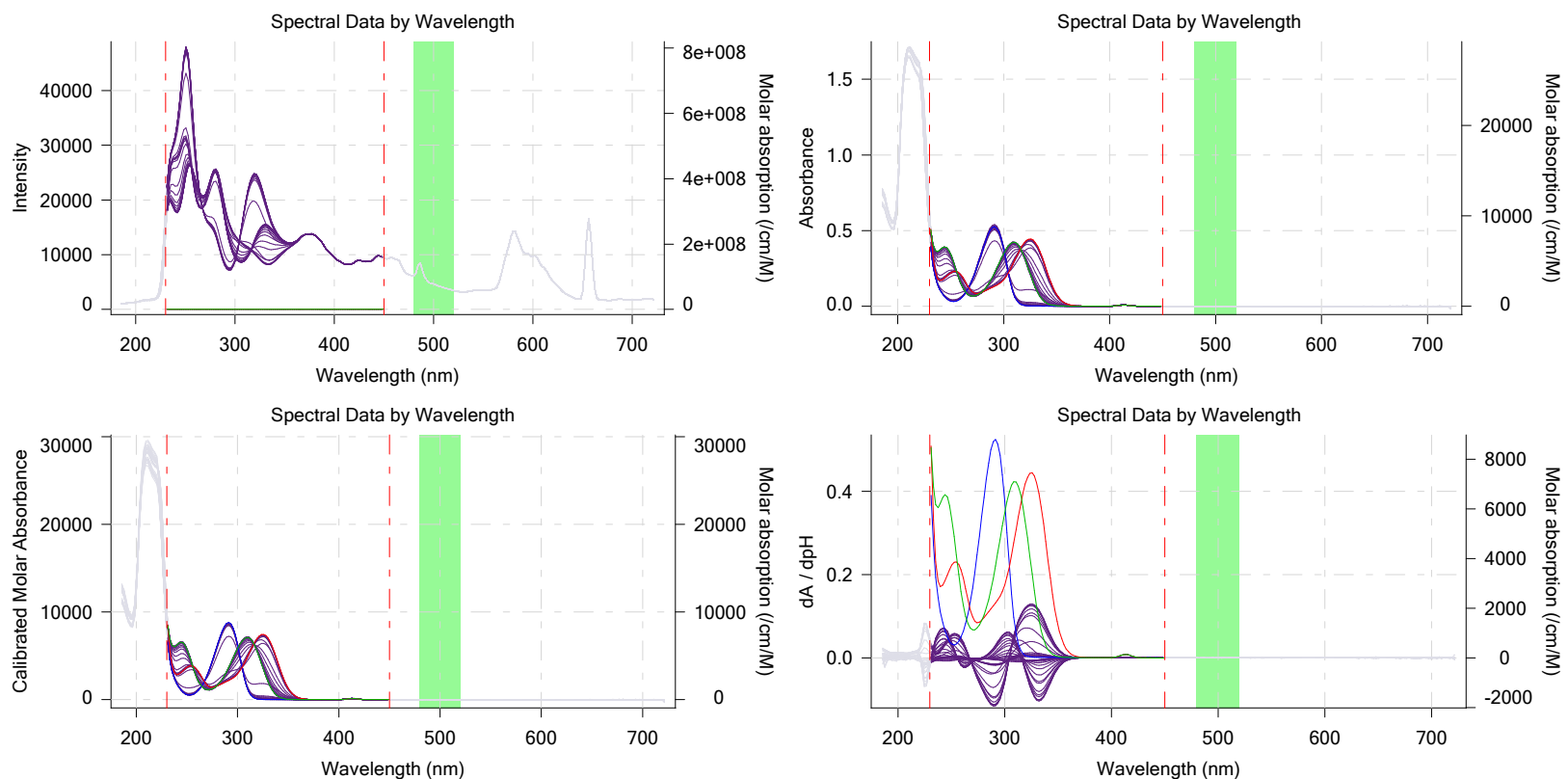

Sample name: **Pyridoxine HCl**  
 Assay name: **UV-metric pKa**  
 Assay ID: **18E-22004**  
 Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

Experiment start time: **5/22/2018 10:11:10 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**


## Graphs (continued)

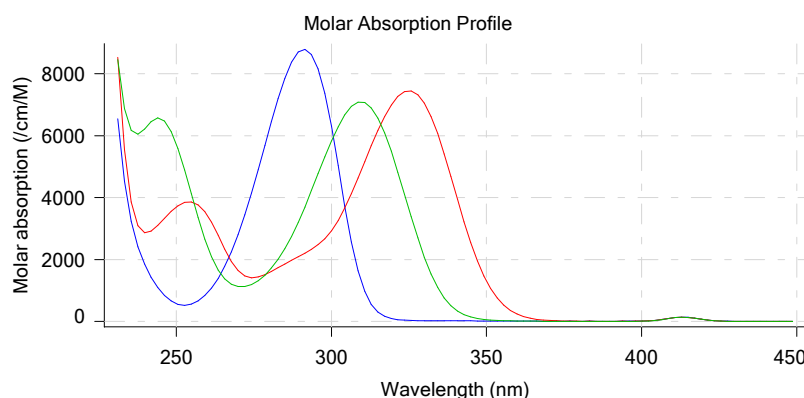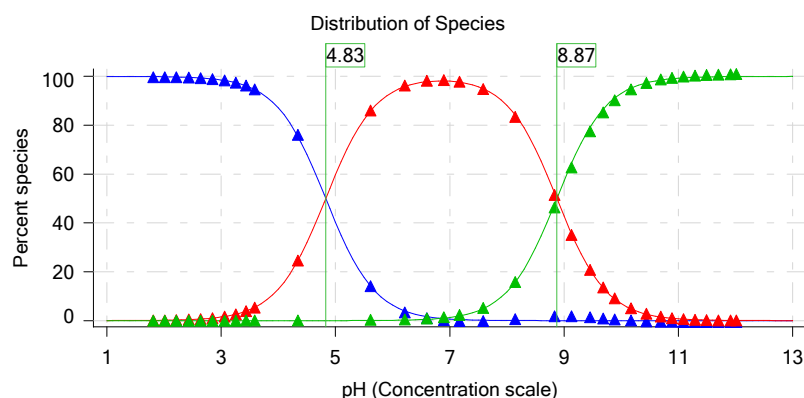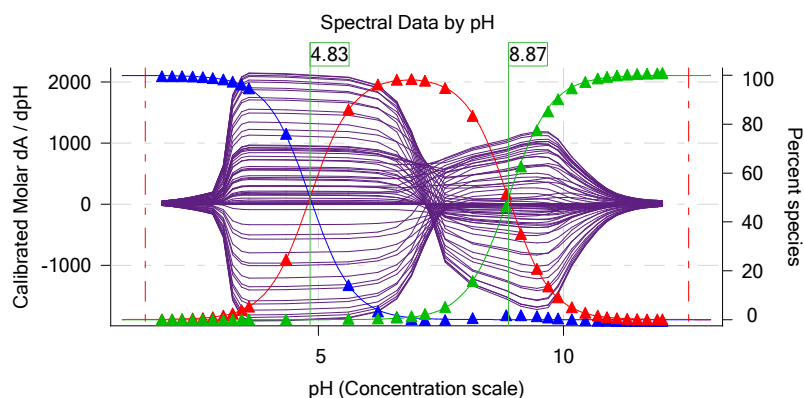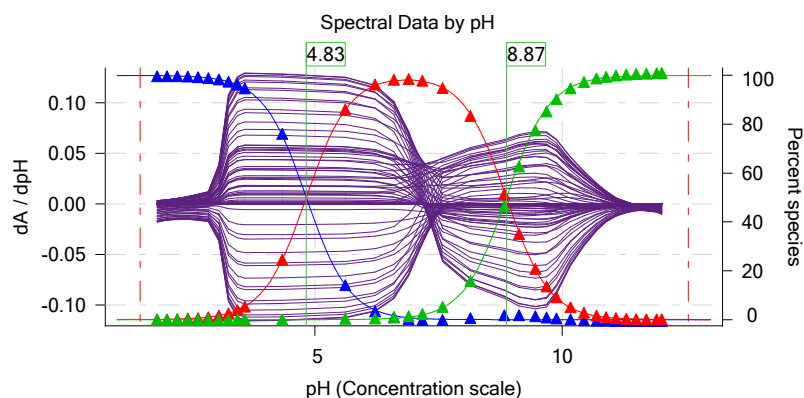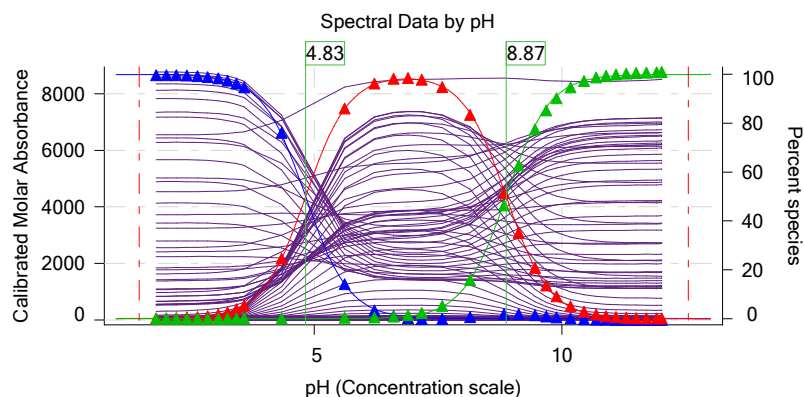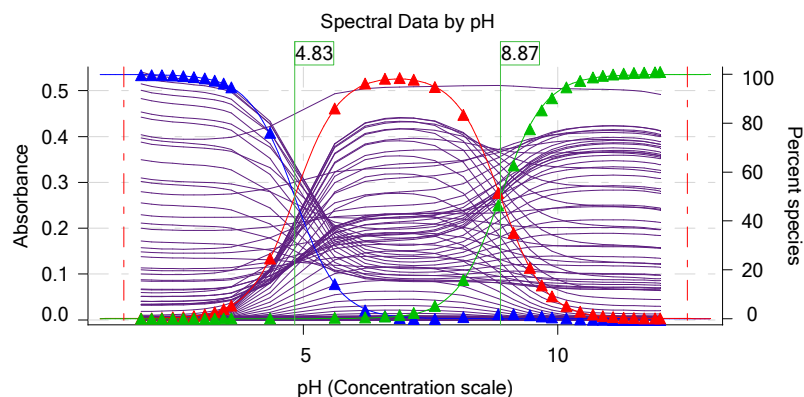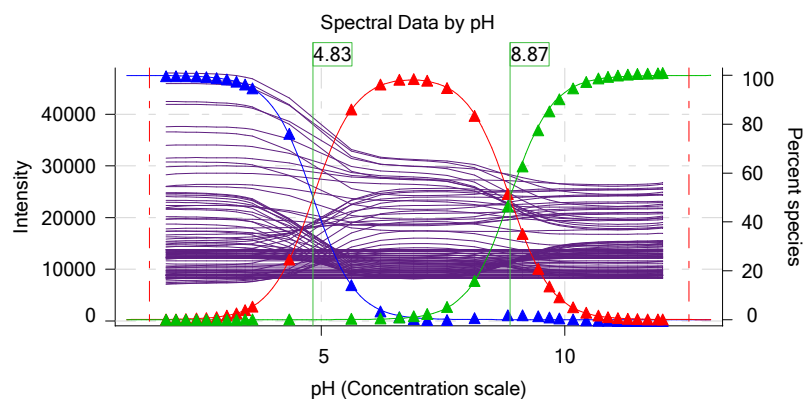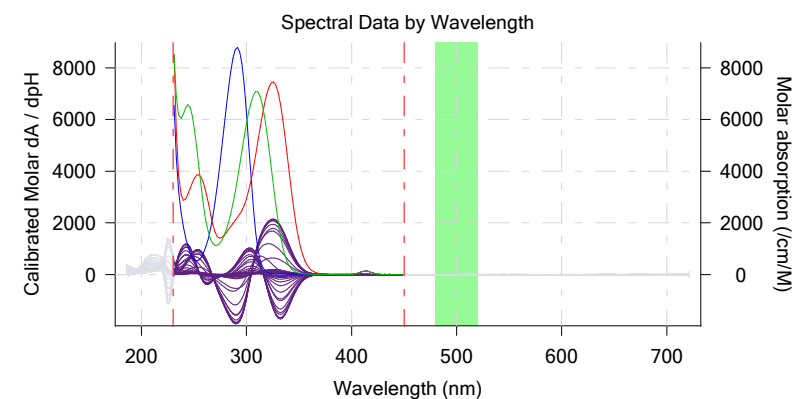

Sample name: **Pyridoxine HCl**  
Assay name: **UV-metric pKa**  
Assay ID: **18E-22004**  
Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

Experiment start time: **5/22/2018 10:11:10 AM**  
Analyst: **Dorothy Levorse**  
Instrument ID: **T311053**

## Graphs (continued)

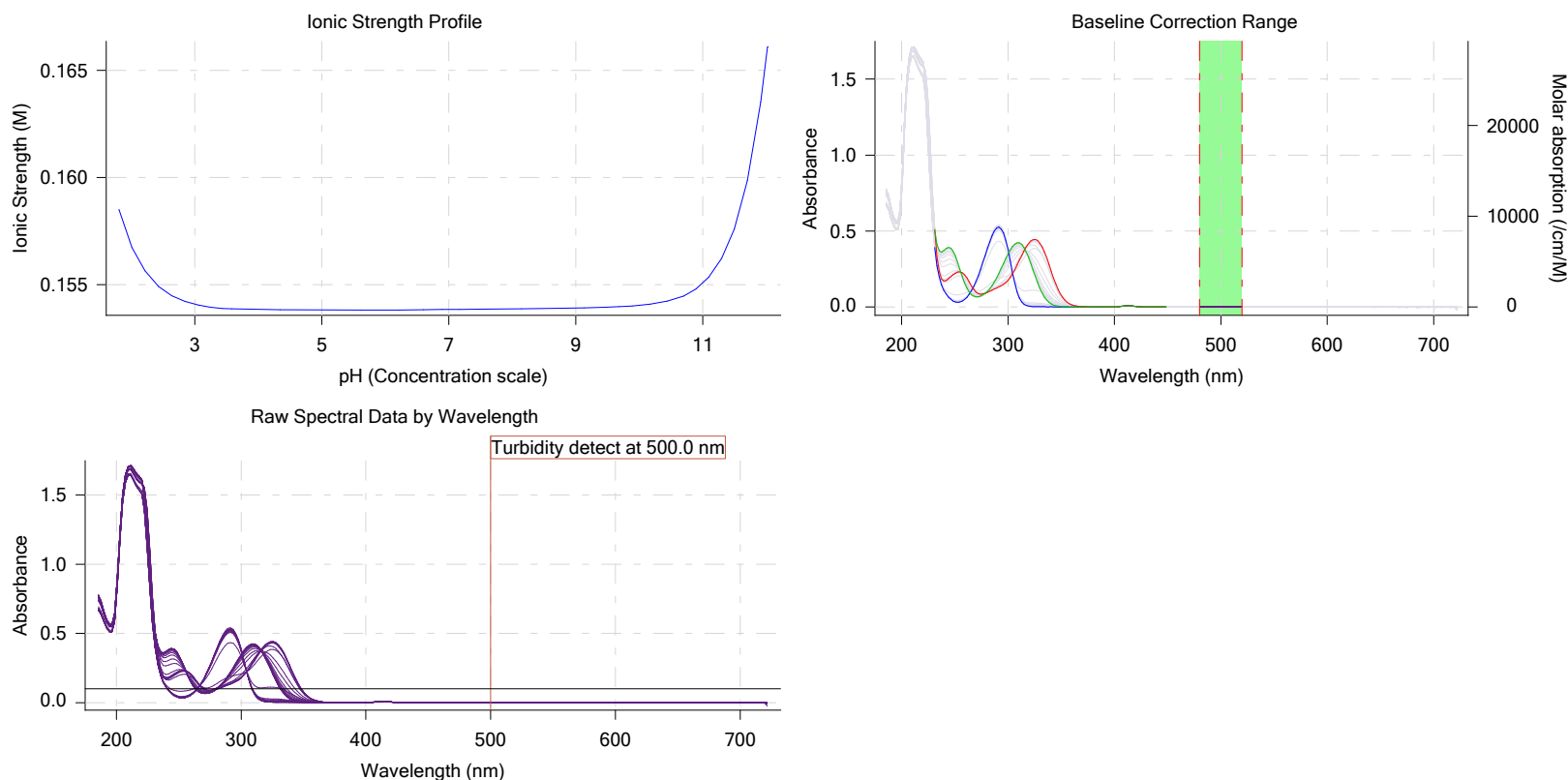

## Assay Model

Settings	Value	Date/Time changed	Imported from
Sample name	Pyridoxine HCl	5/22/2018 9:07:27 AM	User entered value
Sample by	Volume		Default value
Sample volume	0.0020 mL	5/22/2018 9:07:27 AM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.048630 M	5/22/2018 9:07:27 AM	User entered value
Solubility	Unknown		Default value
Molecular weight	205.64	5/22/2018 9:07:35 AM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	2	5/22/2018 9:07:27 AM	User entered value
Sample is a	Ampholyte	5/22/2018 9:07:27 AM	User entered value
pKa 1	4.90	5/22/2018 9:07:27 AM	User entered value
Type	Base	5/22/2018 9:07:27 AM	User entered value
pKa 2	8.80	5/22/2018 9:07:27 AM	User entered value
Type	Acid	5/22/2018 9:07:27 AM	User entered value
logp (XH2 +)	-10.00		Default value
logP (neutral XH)	-10.00	5/22/2018 9:07:27 AM	User entered value
logP (X -)	-10.00		Default value
Stoichiometry	1.00000		Default value
Aprotic counterion name	Chloride		From standards.xml file
Stoichiometry	1.00		From standards.xml file
Charge per counterion	-1		From standards.xml file

## Events

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD	dpH/dt time	Temperature	Intensity Deviation
3:24.7	Dark spectrum											
3:26.3	Reference spectrum											

Sample name: **Pyridoxine HCl**  
 Assay name: **UV-metric pKa**  
 Assay ID: **18E-22004**  
 Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

Experiment start time: **5/22/2018 10:11:10 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Events (continued)

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD
3:54.0	Volume reset due to vial change								
5:24.4	Initial pH = 7.03								
6:27.6	Data point 4	1.50000 mL	0.04986 mL	0.00000 mL	0.02500 mL	1.966	-0.00757	0.56671	0.000
6:56.3	Data point 5	1.50000 mL	0.04986 mL	0.01804 mL	0.02500 mL	2.164	0.00943	0.62450	0.000
7:13.2	Data point 6	1.50000 mL	0.04986 mL	0.02944 mL	0.02500 mL	2.360	-0.00099	0.01699	0.000
7:30.1	Data point 7	1.50000 mL	0.04986 mL	0.03669 mL	0.02500 mL	2.571	-0.00281	0.05908	0.000
7:46.9	Data point 8	1.50000 mL	0.04986 mL	0.04111 mL	0.02500 mL	2.777	0.01674	0.88270	0.000
8:03.6	Data point 9	1.50000 mL	0.04986 mL	0.04381 mL	0.02500 mL	2.981	0.02333	0.92685	0.001
8:20.3	Data point 10	1.50000 mL	0.04986 mL	0.04548 mL	0.02500 mL	3.192	0.01561	0.92316	0.000
8:37.0	Data point 11	1.50000 mL	0.04986 mL	0.04652 mL	0.02500 mL	3.382	0.01511	0.89366	0.000
8:53.7	Data point 12	1.50000 mL	0.04986 mL	0.04718 mL	0.02500 mL	3.561	0.02978	0.98482	0.001
9:10.4	Data point 13	1.50000 mL	0.04986 mL	0.04760 mL	0.02500 mL	3.710	0.05023	0.98647	0.002
9:37.4	Data point 14	1.50000 mL	0.04986 mL	0.04873 mL	0.02500 mL	4.460	0.09744	0.98898	0.004
10:26.4	Data point 15	1.50000 mL	0.04986 mL	0.04915 mL	0.02500 mL	5.719	0.09590	0.94209	0.004
11:29.9	Data point 16	1.50000 mL	0.04986 mL	0.04944 mL	0.02500 mL	6.315	0.09835	0.99059	0.004
12:05.7	Data point 17	1.50000 mL	0.04986 mL	0.04965 mL	0.02500 mL	6.697	0.09733	0.95539	0.004
12:40.6	Data point 18	1.50000 mL	0.04986 mL	0.04981 mL	0.02500 mL	6.990	0.09845	0.98806	0.004
13:14.5	Data point 19	1.50000 mL	0.04986 mL	0.04993 mL	0.02500 mL	7.264	0.09988	0.99331	0.004
13:59.0	Data point 20	1.50000 mL	0.04986 mL	0.05005 mL	0.02500 mL	7.671	0.09246	0.95444	0.004
14:57.4	Data point 21	1.50000 mL	0.04986 mL	0.05016 mL	0.02500 mL	8.226	0.09551	0.91125	0.004
15:49.3	Data point 22	1.50000 mL	0.04986 mL	0.05033 mL	0.02500 mL	8.907	0.09884	0.99191	0.004
16:29.6	Data point 23	1.50000 mL	0.04986 mL	0.05045 mL	0.02500 mL	9.204	0.10049	0.98916	0.004
16:59.7	Data point 24	1.50000 mL	0.04986 mL	0.05061 mL	0.02500 mL	9.525	0.09602	0.98217	0.004
17:24.8	Data point 25	1.50000 mL	0.04986 mL	0.05078 mL	0.02500 mL	9.751	0.09352	0.96911	0.004
17:47.1	Data point 26	1.50000 mL	0.04986 mL	0.05096 mL	0.02500 mL	9.953	0.07957	0.96208	0.004
18:03.9	Data point 27	1.50000 mL	0.04986 mL	0.05125 mL	0.02500 mL	10.234	0.01091	0.87155	0.000
18:30.9	Data point 28	1.50000 mL	0.04986 mL	0.05181 mL	0.02500 mL	10.504	-0.00487	0.70957	0.000
18:47.7	Data point 29	1.50000 mL	0.04986 mL	0.05273 mL	0.02500 mL	10.751	0.00034	0.00611	0.000
19:15.8	Data point 30	1.50000 mL	0.04986 mL	0.05430 mL	0.02500 mL	10.945	-0.00631	0.51904	0.000
19:32.5	Data point 31	1.50000 mL	0.04986 mL	0.05684 mL	0.02500 mL	11.143	0.00276	0.28131	0.000
19:49.4	Data point 32	1.50000 mL	0.04986 mL	0.06087 mL	0.02500 mL	11.345	-0.00336	0.25554	0.000
20:06.3	Data point 33	1.50000 mL	0.04986 mL	0.06738 mL	0.02500 mL	11.544	0.00542	0.72199	0.000
20:23.3	Data point 34	1.50000 mL	0.04986 mL	0.07780 mL	0.02500 mL	11.747	0.00335	0.20625	0.000
20:40.4	Data point 35	1.50000 mL	0.04986 mL	0.09473 mL	0.02500 mL	11.951	0.01060	0.78631	0.000
20:57.3	Data point 36	1.50000 mL	0.04986 mL	0.10687 mL	0.02500 mL	12.052	0.00724	0.77837	0.000
22:56.3	Assay volumes	1.75000 mL	0.16853 mL	0.10687 mL	0.02500 mL				

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
<b>General Settings</b>				
Analyst name	Dorothy Levorse			
Separate reference vial	Yes			
<b>Standard Experiment Settings</b>				
Number of titrations	1			
Minimum pH	2.000			
Maximum pH	12.000			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			
Start titration using	Cautious pH adjust			
<b>Advanced General Settings</b>				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			



## Assay Settings

Sample name: **Pyridoxine HCl**  
Assay name: **UV-metric pKa**  
Assay ID: **18E-22004**  
Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

Experiment start time: **5/22/2018 10:11:10 AM**  
Analyst: **Dorothy Levorse**  
Instrument ID: **T311053**

### Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
<b><i>Titrant Pre-Dose</i></b>				
Titrant pre-dose	None			
<b><i>Assay Medium</i></b>				
Cosolvent in use	No			
ISA water volume	1.50 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
<b><i>Sample Sonication</i></b>				
Sonicate	No			
<b><i>Sample Dissolution</i></b>				
Perform a dissolution stage	No			
<b><i>Carbonate purge</i></b>				
Perform a carbonate purge	No			
<b><i>Temperature Control</i></b>				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
<b><i>Titration 1</i></b>				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
<b><i>Data Point Stability</i></b>				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
<b><i>Experiment cleanup</i></b>				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

### Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.160	5/22/2018 10:11:10 AM	C:\Sirius_T3\18E-21006_Blank standardisation.t3r
Four-Plus S	0.9904	5/22/2018 10:11:10 AM	C:\Sirius_T3\18E-21006_Blank standardisation.t3r
Four-Plus jH	1.0	5/22/2018 10:11:10 AM	C:\Sirius_T3\18E-21006_Blank standardisation.t3r
Four-Plus jOH	-0.5	5/22/2018 10:11:10 AM	C:\Sirius_T3\18E-21006_Blank standardisation.t3r
Base concentration factor	1.012	5/22/2018 10:11:10 AM	C:\Sirius_T3\KOH18D10.t3r
Acid concentration factor	1.011	5/22/2018 10:11:10 AM	C:\Sirius_T3\18E-21006_Blank standardisation.t3r



Sample name: **Pyridoxine HCl**  
 Assay name: **UV-metric pKa**  
 Assay ID: **18E-22004**  
 Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

Experiment start time: **5/22/2018 10:11:10 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 6:24:52 AM
Dispenser 0	Water		3/31/2009 6:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	2-6-18	5/15/2018 2:12:22 PM
Dispenser 2	Acid		3/31/2009 6:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	3-22-18	5/15/2018 2:12:48 PM
Dispenser 1	Base		3/31/2009 6:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	3-22-18	5/15/2018 2:12:34 PM
Dispenser 5	Cosolvent		3/31/2009 6:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 6:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	2-8-18	5/15/2018 2:14:14 PM
Port B	Cyclohexane		4/10/2018 8:40:51 AM
Port C	MeCN (50%, 0.15 M KCl)	4-16-18	5/15/2018 2:14:20 PM
Dispenser 3	Buffer		8/3/2010 6:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Dodecane	1-31-2018	5/15/2018 2:12:54 PM
Dispenser 6	Octanol		10/22/2010 11:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	1-31-2018	4/9/2018 9:14:11 AM
Titration		T3TM1100153	3/31/2009 6:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 10:21:54 AM
E0 calibration	-9.08 mV		5/22/2018 10:11:34 AM
Filling solution	3M KCl	KCL095	5/21/2018 8:57:01 AM
Liquids			
Wash 1	50% IPA:50% Water		5/22/2018 8:38:15 AM
Wash 2	0.5% Triton X-100 in H2O		5/22/2018 8:38:18 AM
Buffer position 1	pH7 Wash		5/22/2018 8:38:22 AM
Buffer position 2	pH 7		5/22/2018 8:38:25 AM
Storage position			5/22/2018 8:38:32 AM
Wash water	3.8e+003 mL	5-15-18	5/15/2018 2:11:48 PM
Waste	6.7e+003 mL		3/19/2018 10:48:12 AM
Temperature controller			8/5/2010 7:35:13 AM
Turbidity detector			3/31/2009 6:24:45 AM
Spectrometer		072390	11/23/2010 12:22:28 PM
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	897:26:49		11/23/2010 12:22:28 PM
Calibrated on	5/21/2018 2:44:22 PM		



Sample name: **Pyridoxine HCl** Experiment start time: **5/22/2018 10:11:10 AM**  
Assay name: **UV-metric pKa** Analyst: **Dorothy Levorse**  
Assay ID: **18E-22004** Instrument ID: **T311053**  
Filename: **C:\Sirius\_T3\18E-22004\_Pyridoxine HCl\_UV-metric pKa\_0417936-0002.t3r**

**Instrument Settings (continued)**

Setting	Value	Batch Id	Install date
Integration time	19		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 10:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

**Refinement Settings**

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050